



TABLE D'OPERATION ELECTRIQUE

Table d'opération TRD-01


The electrically operated surgery table for execution common surgeries and so on chest, abdomen surgical, ophthalmology, otolaryngology, gynecology and obstetrics, urology uses. This product has the unique merit:

The floor may make the longitudinal shift, the Taiwan kneading board uses the high strength plastic film to make, may carry on the beam with the C arm necessary use diagnostic or the shooting;

The leg board is dismountable, manual revolving abducts, under folds, the adjustment is convenient, very convenient urology department surgery;

Grasps the manipulator to use the 24V DC voltage easily, the operation, safe reliable;

Before the floor fluctuation, the caster, about leans and so on main body posture adjustment by the pressed key operation, the electrically operated throughout lever transmission to realize.

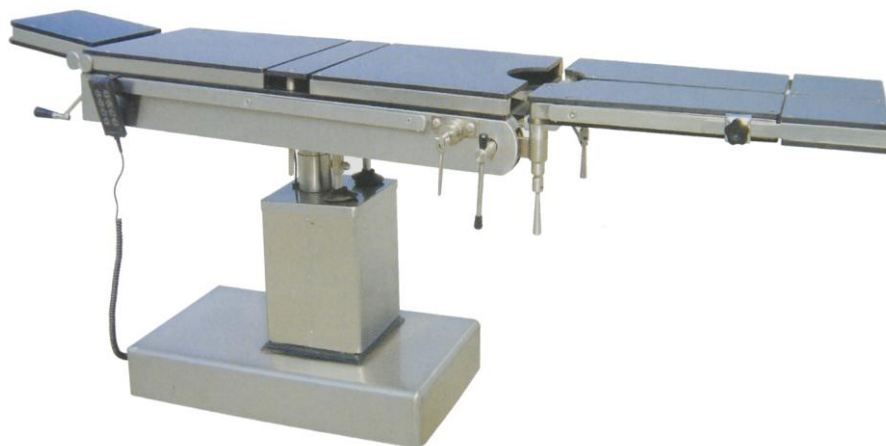

SPECIFICATIONS:

Dimensions of Table Top:

Overall length.....	2100mm
Width.....	480mm

Dimensions of sections:

head board.....	280x300mm
back section.....	500x480mm
seat section.....	450x480mm
leg board.....	500x230mm
foot plate... ..	200x260mm
waist board.....	100x480mm

Table Top from Floor:

Max. Height.....	≤1000mm
Min. Height.....	≥750mm

Adjustment Scope of the Table Top:

forward and backward motion.....	≥25°
leaning from left to right.....	≥20°
longitudinal motion.....	≥300mm

Adjustment Scope of Sections:

Head section raised from the horizon.....	≥45°
Head section lowered from the horizon.....	≥90°
Ascending of the waist board.....	≥100mm
Leg board lowered from the horizon.....	≥90°
Outreaching of the leg board.....	≥90°
The foot boards can be in horizontal level or be vertical with the leg boards	