



UNITE ORL MINI 710

MINI E.N.T TREATMENT UNIT

Twin Oil-less Pump Design
Vacuum Pump & Pressure Pump



CE 0434
ISO 9001
ISO 13485

Designed For Better Life



UNITE ORL MINI 710

Explanation of control panel



- 1 Spray regulator
- 2 Suction regulator
- 3 Suction light
- 4 Spray light
- 5 Lamp switch
- 6 Power switch
- 7 Vacuum gauge
- 8 Spray gauge
- 9 Suction connector
- 10 Spray gun 3 sets

Explanation of position of units




- A Stainless steel round large cotton jar is used to store gauze, cotton balls, etc.
- B Stainless steel square tray with cover is used to store nose mirrors, metal tongue depressors, throat cotton rods, ear cotton rods, etc.
- C The tongue depressor has 3 storage slots for storing tongue depressor, ear cotton rod, etc.
- D Four spare liquid medicine bottles (30cc) are used to store spare liquid medicine bottle, medicine, etc.

Explanation of two sides of instrument



- One 3000cc suction bottle at the right side of the instrument is equipped with automatic full water stop control to prevent dirty water flowing into the pump.
- The suction bottle is made of plastic, can resist hit and high temperature & pressure disinfection. But any chemicals such as antiseptic solution, cleaner, etc. can't be put in the bottle to avoid bottle breaking.
- Installed with two plug sockets for conveniently using together with other instrument.
- PU anti-static moving wheels features with diagonal brake & fixing function.

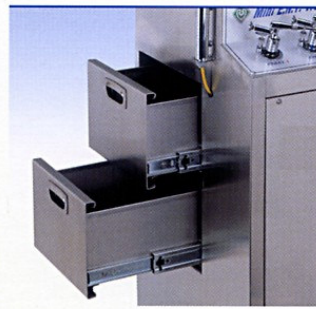

- Two article drawers (L35 x W20 x H12cm) at the left side of the instrument can be used to store articles medicine, catheters, etc.
- The drawer adopts the sliding rail design and can be fully opened for easily taking articles.

STANDARD ACCESSORIES			
Metal suction head		3pcs	
Nose mirrors (big/small)		1pc	
Throat cotton rod		10pcs	
Fuse		2pcs	
SPECIFICATIONS			
Voltage	AC 120V or 240V	Compressor Pump	200W
Current	5A	Max Pressure	1.0~1.3 kg/cm ²
Bulb No.	1	Spray Gun	3
Light Bulb	12V 35W	N.W.	50 KGS
Vacuum Pump	200W	G.W.	60 KGS
Max suction	680~720 mm/hg	Packing Size	L57 X W41 X H92 cm
Max Volume	3,000c.c.		