



AUTOCLAVE 140 à 340 LITRES



AUTOMATIC AUTOCLAVE

Horizontal / Free Stand Type

*Automatic Cycle Control by Electric Process Timers & Solenoid Valves
from Heating Up » Sterilizing » Drying » Exhaust » Cycle End*


Exist in 3 models: 141 - 196 - 339 LITERS




CE 0434
ISO 9001
ISO 13485



AUTOCLAVE 140 à 340 LITRES

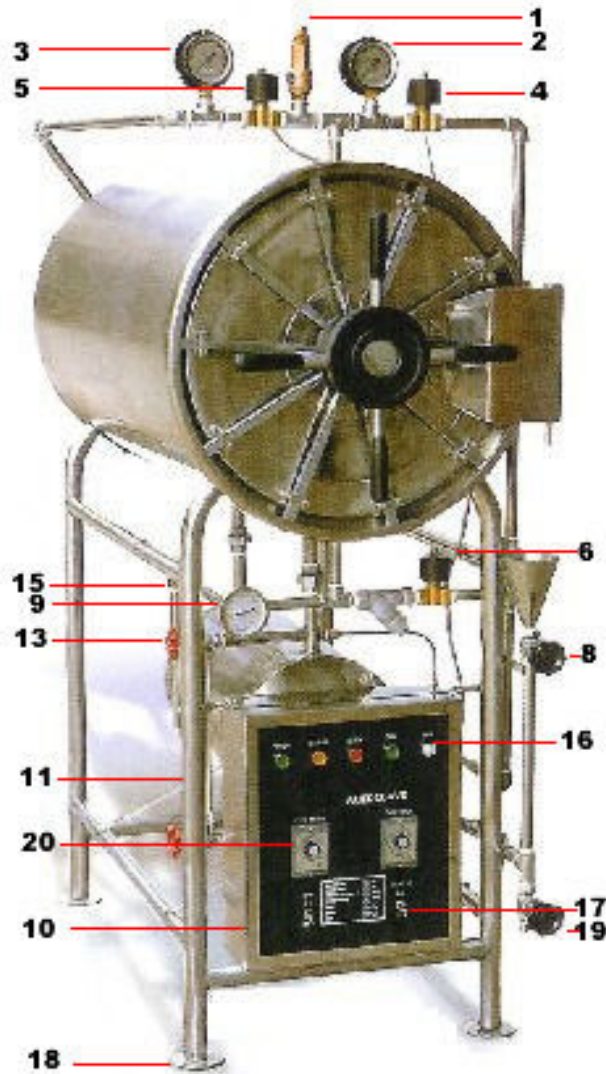

Automatic cycle control by electric process timers & solenoid valves.

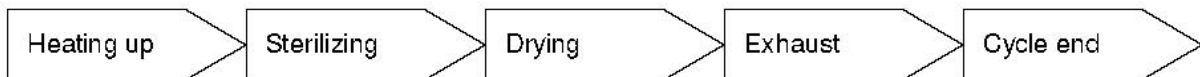



AUTOCLAVE 140 à 340 LITRES

INSTRUMENTS & MOUNTINGS			
1	Safety release valve	11	Water level gauge
2	Jacket pressure gauge	12	Automatic pressure controller
3	Chamber pressure gauge	13	Steam generator
4	Drying valve	14	Automatic "Burn-out" protection
5	Sterilizing valve	15	Steam trap
6	Exhaust valve	16	Heater pilot lamp
7	Chamber bottom plate. (Stainless steel)	17	Main switch ON/OFF
8	Water fill valve	18	Level adjustable flange
9	Dial thermometer	19	Water drain valve
10	Electrical control box	20	Adjustable timer
STANDARD ACCESSORIES			
Heating element		1pc	
Silicone door seal		1pc	

SPECIFICATIONS MODEL 141 L			
Power Supply	AC 220/380V, 7 KW	Chamber diameter	45 cm
Sterilizing pressure	1.2~2.1 kg/cm \diamond	Chamber depth	90 cm
Sterilizing temperature	121~134 _ C	N.W. / B.W.	210 KGS / 300 KGS
Chamber capacity	141 liters	Packing size	L137 x W79 x H187 cm

SPECIFICATIONS MODEL 196 L			
Power Supply	AC 220/380V, 9 KW	Chamber diameter	50 cm
Sterilizing pressure	1.2~2.1 kg/cm \diamond	Chamber depth	100 cm
Sterilizing temperature	121~134°C	N.W. / B.W.	240 KGS / 340 KGS
Chamber capacity	196 liters	Packing size	L148 x W85 x H187 cm

SPECIFICATIONS MODEL 339 L			
Power Supply	AC 220/380V, 12 KW	Chamber diameter	60 cm
Sterilizing pressure	1.2~2.1 kg/cm \diamond	Chamber depth	120 cm
Sterilizing temperature	121~134 _ C	N.W. / B.W.	350 KGS / 400 KGS
Chamber capacity	339 liters	Packing size	L167 x W96 x H193 cm